



**MILESTONE**  
H E L P I N G  
P A T I E N T S

# MileWATCH

## 24/7 Unit surveillance system



MileWATCH is a 24/7 units surveillance appliance that continuously monitors and controls Milestone units throughout day and night. Status Information are collected from units installed in both Operating Theaters and Pathology Laboratory. This, in turn, also helps having a better planning of the Laboratory daily work. Compatible with existing Laboratory PC and Server, MileWATCH works without the need of an internet connection: all the related information is securely stored within the local LAN.


MileWATCH offers the possibility to create a unique samples' Chain of Custody, from the Operation Theater, preserving evidence of the case ID, formalin ratio used, temperature and time of the transportation, until tissues have arrived to the Pathology Lab. Tissue processing protocol details and final storage information info will also be collected by MileWATCH and sent to local LIS (Laboratory Information System).



- » **MONITOR THE PROCESS IN REAL TIME**
- » **BE NOTIFIED BY E-MAIL**
- » **GET EVENTS LOG AND LOG FILE**
- » **SORT DOCUMENTATION BY PATIENT-CASE ID**
- » **INTEGRATE INTO THE LIS**

# Remote monitoring

MileWATCH software provides extended instrument control from outside the laboratory with live view of the units screen. This allows the remote monitoring of every process status during night and day. An instant audio visual alert provides assistance 24 hours a day, 7 days a week. Moreover, all units are connected to server for a prompt transfer of process documentation. All documentation is related to the hospital Patient-Case ID.

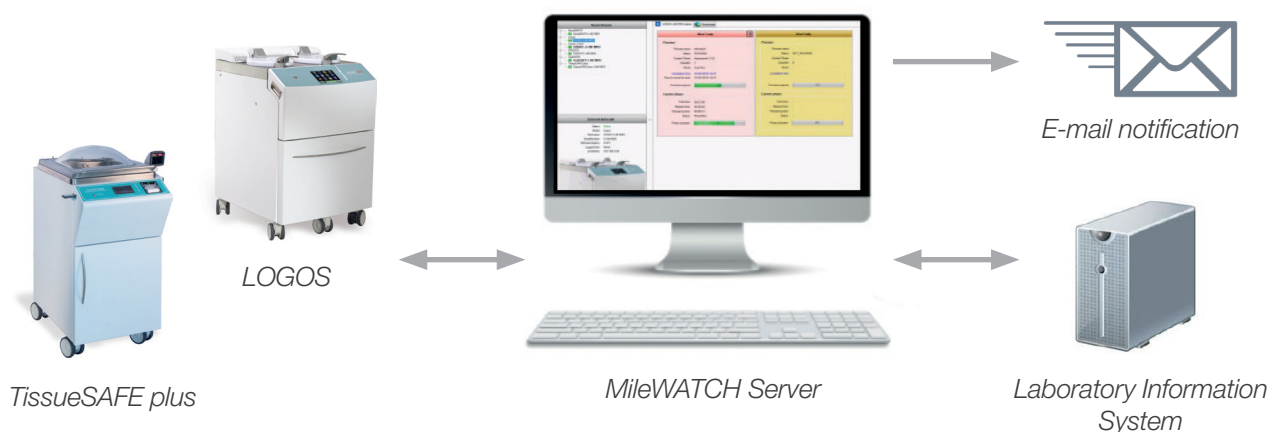

## E-mail notification system

MileWATCH can be set to send live e-mail notification to the interested parties providing all the necessary information for a prompt action to assure trouble-free operations. Recipients can be customized based on notification type. At the end of the program, an e-mail automatically notifies mayor occurrences (events log and log file) to the service personnel and protocol status to the laboratory manager.

## LIS integration

Accurate and precise process documentation can be automatically collected and directly transferred by interfacing MileWATCH to the LIS (Laboratory Information System) through the local LAN. It can be stored and retrieved, sorting by Case ID.

The incredibly intuitive MileWATCH software further enhances the innovative features of the Milestone units, and makes this tool the best choice to answer to the need for documentation and tracking.

*MileWATCH is available for LOGOS series, TissueSAFE plus, and SealSAFE, Temperature Data Logger.*

### Headquarters

MILESTONE Srl  
Via Fatebenefratelli, 1/5  
24010 Sorisole (BG) - Italy  
Tel: +39 035 4128264  
Fax: +39 035 575498  
email: [marketing@milestonemed.com](mailto:marketing@milestonemed.com)

### US Office

MILESTONE MEDICAL TECHNOLOGIES, INC.  
6475 Technology Avenue  
Suite F, Kalamazoo, MI 49009 - USA  
Tel: 269-488-4950 - Toll-free: 866-995-5300  
Fax: 269-488-4949  
email: [info@milestonemed.com](mailto:info@milestonemed.com)