



MILESTONE
H E L P I N G
P A T I E N T S



WorkSTATION


**Mobile Grossing Station
for Small Biopsies**



Small Biopsy Grossing

WorkSTATION

WorkSTATION is an all-in-one grossing station for small biopsies. Ergonomically designed to standardize grossing procedures for a more efficient and complete diagnostic process, WorkSTATION ensures improved safety for both the operator and specimen. WorkSTATION includes features such as easily adjustable working height, an all-in-one control panel, and backdraft and downdraft exhaust. Paired with Milestone's MacroPATH digital imaging system, WorkSTATION provides complete and standardized digital documentation, including on-screen annotations and tags, high quality images, and options for voice and video documentation.

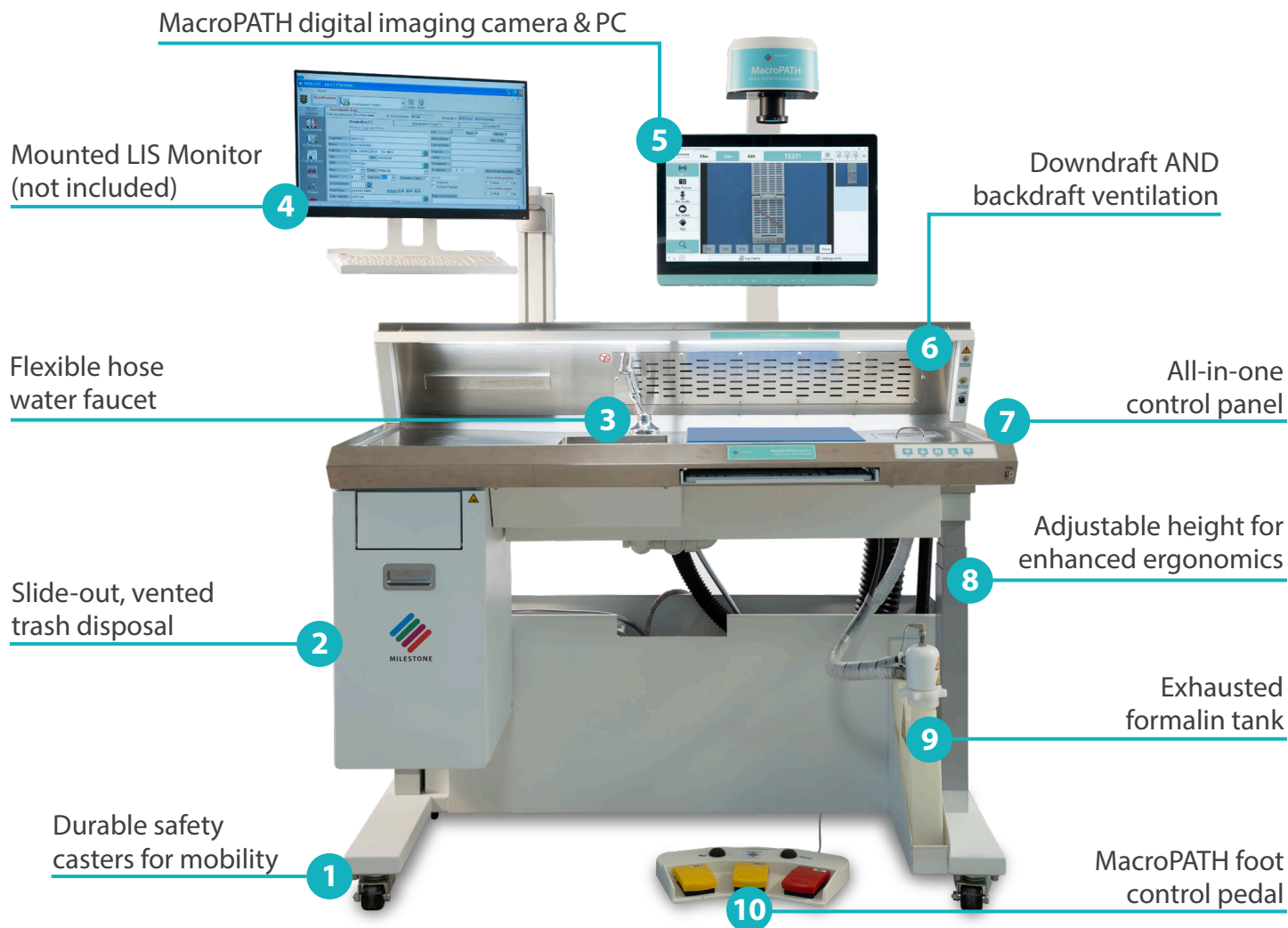

Milestone's Line of Grossing Stations

WorkSTATION



eGROSS



UltraGROSS



THE POWER OF DIGITAL IMAGING

WorkSTATION and Milestone's MacroPATH digital imaging system are a perfect pair. With MacroPATH, Pathologists' Assistants and grossing technicians can digitally document the grossing process by taking images, videos, and voice recordings to provide Pathologists with critical data for evaluation and final diagnosis.



- ✓ High Resolution Camera
For recording images and videos
- ✓ Touch-Screen Monitor
For easy and convenient operation
- ✓ Tagging Function
For easy case/file retrieval
- ✓ Remote Viewing
For easy Pathology consultations
- ✓ Suite of Editing Tools
For on-screen sizing and editing

KEY BENEFITS

OPTIMAL ERGONOMICS

- Optimal working position
- Motorized, height adjustable working area
- Casters with safety locks for mobility

ENHANCED USER SAFETY

- Powerful downdraft and backdraft exhaust
- Vented wastepaper basket
- Designed to minimize formalin exposure

ADVANCED FEATURES

- High resistance durable working surface
- Easily accessible, user friendly control panel
- Foot pedal control for MacroPATH functions

DIGITAL DOCUMENTATION

- Digital documentation with MacroPATH
- Ensures complete traceability/chain of custody
- Full suite of editing tools for on-screen annotations

TRACKING & MONITORING

- Complete chain of custody documentation
- Automatic photo capturing with C-Scan (optional)
- Integration of MacroPATH with your LIS

TEACHING TOOL

- Create a database of cases including photos/videos
- Enhance quality of diagnostic reporting
- Full documentation for digital presentations

UNIQUE FEATURES



MACRO DIGITAL DOCUMENTATION

Full documentation is ensured via the built-in MacroPATH macro digital imaging system (optional), the most widely installed imaging system in the world. MacroPATH enhances the quality of diagnostic reports for tumor boards and allows for specimen tracking downstream.



OPTIMAL ERGONOMICS

The working area height can be easily adapted to the position most comfortable for each individual user through the on-board, easy to use, control terminal.



ENHANCED SAFETY

Built for safety with an on-board backdraft and downdraft ventilation system and an exhausted formalin tank with no fume return valve, filter, and light alarm.



EASY ACCESS CONTROL PANEL

Featuring an on-board control panel which provides key functions to optimize the workflow at the touch of a button. Includes height adjustments, exhaust fan and light controls.

RELATED PRODUCTS



MacroPATH:
Digital Imaging
System



MileWATCH:
24/7 Monitoring
& Tracking System



T-Tracker:
Automatic Rack
Scanner



HistoDream TS:
Biospecimen Tracking
System



US Headquarters

MILESTONE MEDICAL TECHNOLOGIES, INC.

6475 Technology Avenue, Suite F - Kalamazoo, MI 49009

Tel: 269-488-4950 - Toll-free: 866-995-5300 - Fax: 269-488-4949

email: info@milestonemed.com - www.milestonemed.com