



**MILESTONE**

H E L P I N G  
P A T I E N T S



# HistoDream EW

Embedding Workstation  
with biospecimens traceability



## HistoDream EW - EMBEDDING WORKSTATION - MAIN FEATURES



---

5 liters of molten paraffin

---

Programmable heating

---

Reservoir & conduits always heated, ready for urgent operations

---

Dispenser compatible with Mega and SuperMega cassette.

---



---

Up to 300 cassette

---

Milestone racks and any standard rack

---

Two built in paraffin trimmers

---



---

Membrane keys

---

Programmable heating including non working holidays

---

Extra touch screen with bar code reader for case traceability

---





---

Capability of up to 400 molds

---

Programmable heating

---

Space suitable for Mega and SuperMega cassette

---

## HistoDream CM - COLD EMBEDDING MODULE - FEATURES:

---

- Large cooling surface
- Selectable temperature
- Steady temperature distribution over the entire surface.
- Stand-alone operation possible.
- Can be placed either on the right or the left side of the dispensing console.



---

Temperature from 0° up to -12° Celsius

---

Space of 60 standard or 45 Synergy molds

---

Environmental safe refrigerant (R134a)

---




## HistoDream EW - EMBEDDING WORKSTATION - INCLUDES:

- HistoDream Dispensing Module
- HistoDream Thermal Module
- HistoDream Cold Embedding Module\*

\*This module can operate also independent

## ENHANCE THE SPECIMENS TRACEABILITY WITH VISUAL

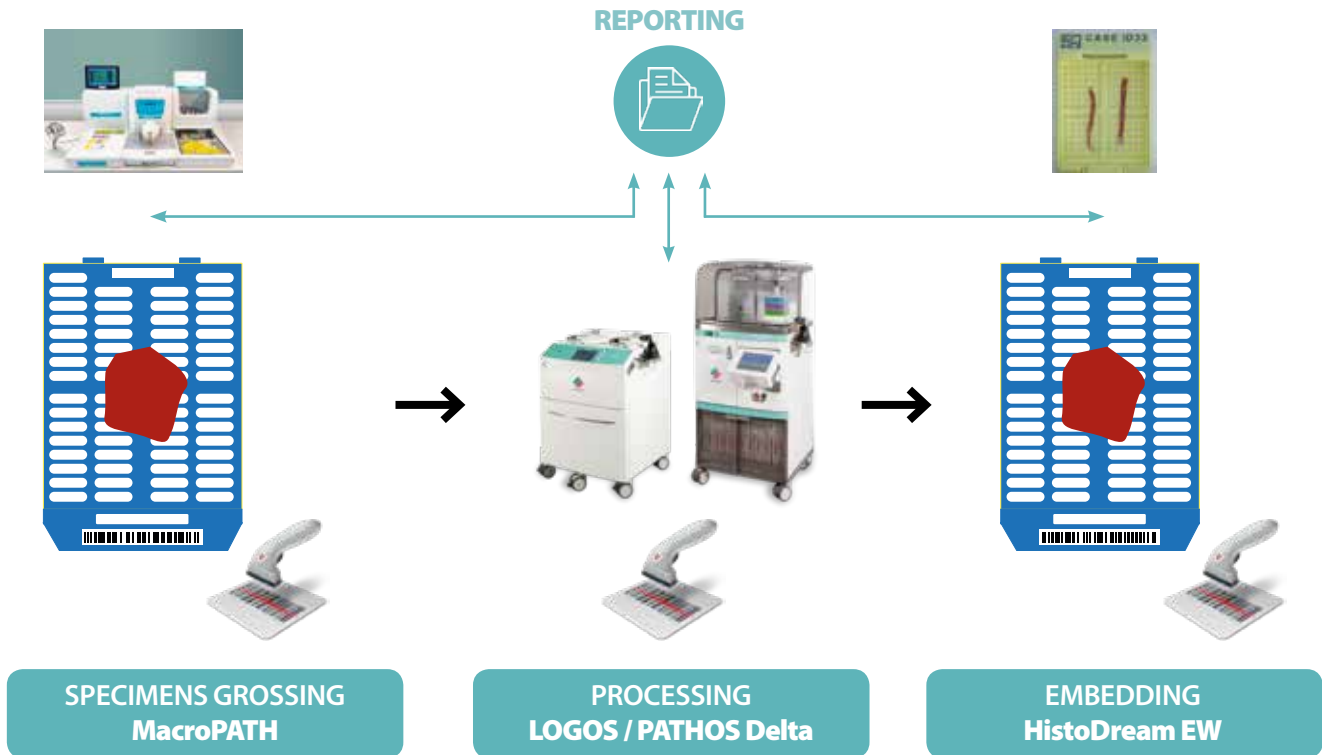

### Optionals

- HistoDream PC: Tracking system terminal, including Viewer Software and 2D bar code reader.
- Magnifying glass.
- Heated tweezers 1 mm, 2mm, 4mm.

### Technical Specifications

- Total Dimensions: 102 cm X 64,1 cm X 40,5 cm (WxDXH).
- Total Weight: 61 Kg.

### Headquarters

MILESTONE Srl  
Via Fatebenefratelli, 1/5  
24010 Sorisole (BG) - Italy  
Tel: +39 035 4128264  
Fax: +39 035 575498  
email: marketing@milestonemed.com

### US Operations

MILESTONE MEDICAL TECHNOLOGIES, INC.  
6475 Technology Avenue  
Suite F, Kalamazoo, MI 49009 - USA  
Tel: 269-488-4950 - Toll-free: 866-995-5300  
Fax: 269-488-4949  
email: info@milestonemed.com